

THE 67<sup>TH</sup> INTERNATIONAL



# OPEN READINGS

CONFERENCE FOR STUDENTS OF PHYSICS AND NATURAL SCIENCES

**BOOK OF  
ABSTRACTS** | **2024**



Vilnius  
University

VILNIUS UNIVERSITY PRESS

Editors:

Martynas Keršys  
Rimantas Naina  
Vincentas Adomaitis  
Emilijus Maskvytis

Cover and Interior Design:

Goda Grybauskaitė

Vilnius University Press  
9 Saulėtekio Av., III Building, LT-10222 Vilnius  
info@leidykla.vu.lt, [www.leidykla.vu.lt/en/](http://www.leidykla.vu.lt/en/)  
[www.knygynas.vu.lt](http://www.knygynas.vu.lt), [www.journals.vu.lt](http://www.journals.vu.lt)

Bibliographic information is available  
on the Lithuanian Integral Library Information System (LIBIS) portal [www.ibiblioteka.lt](http://www.ibiblioteka.lt)  
ISBN 978-609-07-1051-7 (PDF)

© Vilnius University, 2024

## DEER SPECTROSCOPY OF S100A9 PROTEIN

Aistė Peštenytė<sup>1</sup>, Gediminas Usevičius<sup>1</sup>, Darius Šulskis<sup>2</sup>, Ieva Baronaitė<sup>2</sup>, Vytautas Smirnovas<sup>2</sup>, Jūras Banys<sup>1</sup>, Mantas Šimėnas<sup>1</sup>

<sup>1</sup>Faculty of Physics, Vilnius University, Saulėtekio al. 3, LT-10257, Vilnius, Lithuania

<sup>2</sup>Sector of Amyloid Research, Institute of Biotechnology, Life Sciences Centre, Vilnius University, LT-10257 Vilnius, Lithuania  
aiste.pestenyte@ff.stud.vu.lt

Double electron-electron resonance (DEER) spectroscopy, a widely used tool in structural biology, explores biomolecules like proteins, RNA, and DNA by determining nanoscale distances between unpaired electron spins [1]. This pulsed electron paramagnetic resonance (EPR) technique utilizes microwave pulses of varying frequencies. In this process, one electron spin is detected through EPR, while the other one is excited by rotating its magnetization vector (Fig. 1) [2]. Throughout the experiment, we monitor the alterations in dipole interaction between the interacting unpaired electrons. DEER overcomes limitations in studying systems without unpaired electrons, such as biomolecules, by incorporating spin labels through site-directed spin-labeling (SDSL) [2]. Compared to other methods, DEER is not constrained by crystallized molecules or molecular weight limitations. The combination of SDSL with DEER emerges as a highly promising approach for structural analysis, enabling the observation of changes in local structure, interactions with other molecules, alterations in the distance between electron spins, and the visualization of conformational heterogeneity and dynamics [1].

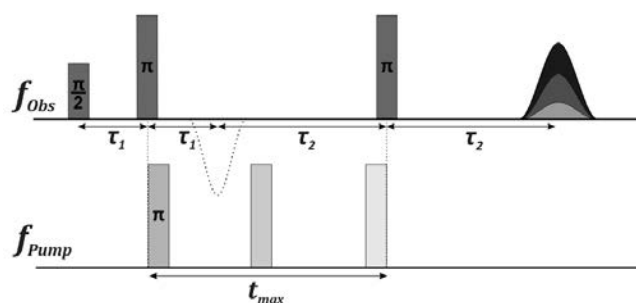


Fig. 1. Microwave pulse sequence of the DEER experiment.

To practically identify capabilities and limitations of DEER spectroscopy, we examined a calcium-binding protein S100A9, related to Alzheimer's and Parkinson's diseases. During this study, we determined the distance distribution (Fig. 2) between the two cysteine groups, utilizing the nitroxide radical as a spin label.

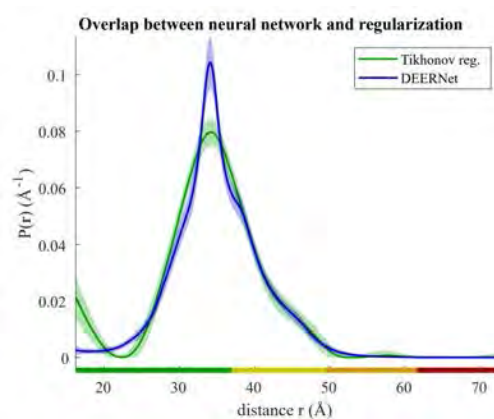


Fig. 2. Distance distribution between electron spins in the S100A9 protein.

[1] Indra D. et al., Use of electron paramagnetic resonance to solve biochemical problems, Biochemistry 2013.

[2] Jeschke G. DEER distance measurements on proteins. Annual review of physical chemistry. 2012 May.